



Model UF V 700 | Ultra-low temperature freezers with climate-neutral refrigerants

The BINDER ultra low temperature freezer ensures the safe storage of samples at -90 °C. It combines outstanding environmental friendliness with low energy consumption, convenient operation and an individual safety concept.

BENEFITS

- · Lowest energy consumption in its class
- Environmentally-friendly refrigerants R-290 and R-170
- Thermal insulation with unique long-life VIPs (vacuum insulation panels)
- · Rust-proof interior made completely of stainless steel



Model UF V 700



Model UF V 700

IMPORTANT FEATURES

- Temperature range: -90 °C to -40 °C
- Lowest energy consumption in its class
- · Powerful cascade compressor cooling unit
- Environmentally friendly refrigerants R-290 and R-170
- Thermal insulation with unique long-life VIPs (vacuum insulation panels)
- Innovative door gasket concept reduces buildup of ice
- Ergonomic door handle
- · Rust-proof interior made completely of stainless steel

- Detachable inner doors made of stainless steel
- 3 stainless-steel shelves providing flexible placement
- Water cooling (depending on version)
- Zero-voltage alarm contact for in-house alarm system
- Internal data logger, measured values can be read out in open format via USB
- · Ethernet interface
- 2 access ports, Ø 28 mm, rear

ORDERING INFORMATION

Designation	Option model	Interior volume	Rated Voltage	Power frequency	Article Number
UFV700-230V	Standard	700 L	230 V	50 Hz	9020-0348
UFV700-230V-W	Water cooling	700 L	230 V	50 Hz	9020-0350
UFV700UL-240V	Standard	700 L	208230 V	6o Hz	9020-0354
UFV700UL-120V	Standard	700 L	115 V	6o Hz	9020-0352




TECHNICAL DATA

Data				
Designation	UFV700-230V	UFV700-230V-W	UFV700UL-240V	UFV700UL-120V
Article Number	9020-0348	9020-0350	9020-0354	9020-0352
Option model	Standard	Water cooling	Standard	Standard
Performance Data Temperature				
Temperature range	-9040 °C	-9040 °C	-9040 °C	-9040 °C
Temperature uniformity at -80°C	2.5 ± K	2.5 ± K	2.5 ± K	2.5 ± K
Temperature fluctuation at -80°C	1.5 ± K	1.5 ± K	1.5 ± K	1.5 ± K
Electrical data				
Rated Voltage	230 V	230 V	208230 V	115 V
Power frequency	50 Hz	50 Hz	60 Hz	60 Hz
Nominal power	1,6 kW	1,6 kW	1,8 kW	1,4 kW
Unit fuse	10,0 A	10,0 A	10,0 A	15,0 A
Phase (Nominal voltage)	1~	1~	2~	1~
Measures				
Interior volume	700 L	700 L	700 L	700 L
Net weight of the unit (empty)	288 kg	288 kg	288 kg	288 kg
Permitted load	200 kg	200 kg	200 kg	200 kg
Permitted load of individual compartment	50 kg	50 kg	50 kg	50 kg
Wall clearance back	100 mm	100 mm	100 mm	100 mm
Wall clearance sidewise	250 mm	250 mm	250 mm	250 mm
Housing dimensions not incl. fittings and connections				
Width net	1,110 mm	1,110 mm	1,110 mm	1,110 mm
Height net	1,966 mm	1,966 mm	1,966 mm	1,966 mm
Depth net	938 mm	938 mm	938 mm	938 mm
Internal Dimensions				
Interior width	890 mm	890 mm	890 mm	890 mm
Interior height	1,300 mm	1,300 mm	1,300 mm	1,300 mm
Interior depth	605 mm	605 mm	605 mm	605 mm
Inner doors	2	2	2	2
Unit doors	1	1	1	1
Environment-specific data				
Sound-pressure level	47 dB(A)	47 dB(A)	47 dB(A)	47 dB(A)
Energieverbrauch bei Sollwert -80 °C und T_u = 20 °C	8.1 kWh/d	8.1 kWh/d	8.1 kWh/d	8.1 kWh/d
Fixtures				
Number of shelves (std./max.)	3/11	3/11	3/11	3/11
Load-variants				
Freezer racks per compartment	6	6	6	6
Cryoboxes, 50 mm	528	528	528	528
All		٠ ، ، ، ، ، ، ، ، ، ، ، ، ، ، ، ، ، ، ،	00 1	I. C:


All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.




DIMENSIONS Incl. fittings and connections [mm]

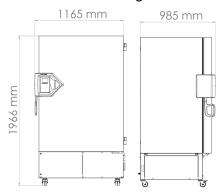

DIAGRAMS

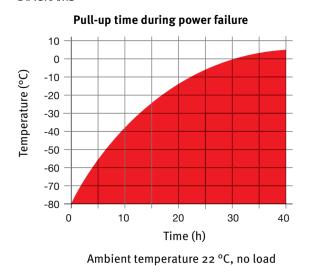

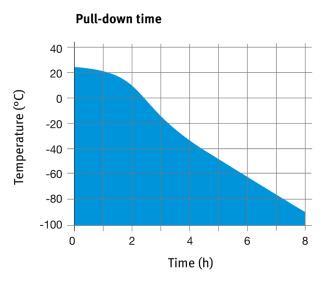

Ambient temperature 25 °C, no load

Heating up rate

Cooling down rate

OPTIONS

Designation	Description	*	Article Number
Alarm system, battery-backed	for maintaining temperature display, alarm messages, and communication		8012-1910
Analog output 4-20 mA	for temperature values (output not adjustable)	_	8012-1898
Buck-boost converter	automatic voltage compensation	_	8012-1897
CO ₂ backup cooling	factory installed (allow for increased depth of UF V: 1085 mm), includes the option of a battery-backed alarm system	_	8012-1936
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	-	8012-1126
Calibration certificate, temperature	for temperature, measurement in center of chamber at specified temperature	_	8012-0031
	temperature measurement incl. certificate, 20 measuring points at specified temperature	⁻ 0	8012-0925
Compartment doors	consisting of electromechanical door locking mechanism and electronic access		8012-1894
Door access system			8012-1899
Door lock	lockable door handle		8012-1911
Shipping pallet	ping pallet with unloading ramp for installing unit without forklift		8012-1907
Switch cover	Lockable main switch cover		8012-1980



ACCESSORIES

Designation	Description	Option model	* Article Number
	C 13 refrigerator connection cable with BS 1363 British plug		5023-0285
5 m power cable	C 13 refrigerator connection cable with CEE 7/7 European plug	_	5023-0283
	C 13 refrigerator connection cable with type 13 Swiss plug	_	5023-0284
APT-COM™ 4 BASIC-Edition	for simple logging and documentation requirements with up to 5 networked units.	version 4, BASIC edition	9053-0039
APT-COM™ 4 GLP-Edition	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.	version 4, GLP edition	9053-0042
APT-COM™ 4 PROFESSIONAL- Edition	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.	version 4, PROFESSIONAL edition	9053-0040
	for sure standard house (so mm) dimensions Home, Wasan Deck, mm	Set of 3	8012-2345
	for 5x4 standard boxes (50 mm), dimensions: H 278 x W 139 x D 564 mm	Set of 4	8012-2346
	for Eyy standard boyos (Fo mm) incl boyos dimensions. H ato y Waso y D 56 / mm	Set of 3	8012-2347
	for 5x4 standard boxes (50 mm) incl. boxes, dimensions: H 278 x W 139 x D 564 mm	Set of 4	8012-2348
Chamber rack, aluminum	for (v. standard bours (s. san) dimensions Heavy West Deck and	Set of 3	8012-2357
	for 6x4 standard boxes (50 mm), dimensions: H 334 x W 139 x D 564 mm	Set of 4	8012-2358
		Set of 3	8012-2359
	for 6x4 standard boxes (50 mm) incl. boxes, dimensions: H 334 x W 139 x D 564 mm	Set of 4	8012-2360
		Set of 3	8012-2349
	for 5x4 standard boxes (50 mm), dimensions: H 278 x W 139 x D 561 mm	Set of 4	8012-2350
		Set of 3	8012-2351
	for 5x4 standard boxes (50 mm) incl. boxes, dimensions: H 278 x W 139 x D 561 mm	Set of 4	8012-2352
Chamber rack, stainless steel		Set of 3	8012-2361
	for 6x4 standard boxes (50 mm), dimensions: H 334 x W 139 x D 561 mm	Set of 4	8012-2362
	for 6x4 standard boxes (50 mm) incl. boxes, dimensions: H 334 x W 139 x D 561 mm	Set of 3	8012-2363
		Set of 4	0 8012-2364
	for 5x4 standard boxes (50 mm), dimensions: H 278 x W 140 x D 562 mm	Set of 3	8012-2353
		Set of 4	8012-2354
	for 5x4 standard boxes (50 mm) incl. boxes, dimensions: H 278 x W 140 x D 562 mm	Set of 3	8012-2355
Chamber rack with drawer,		Set of 4	8012-2356
stainless steel		Set of 3	8012-2365
	for 6x4 standard boxes (50 mm), dimensions: H 334 x W 140 x D 562 mm	Set of 4	8012-2366
		Set of 3	8012-2367
	for 6x4 standard boxes (50 mm) incl. boxes, dimensions: H 334 x W 140 x D 562 mm	Set of 4	
Cryoboxes	set consisting of 36 cardboard boxes with sample dividers	301 01 4	8012-2368 8012-2369
Deicing kit	consisting of 30 cardboard boxes with sample dividers	_	8012-1909
Grid shelf	For temperature range -4o °C to -6o °C, 3 provided in set, including shelf supports	stainless steel	8012-1904
	Q/OQ/PQ documents – supporting documents for validation performed by	Digital in PDF format	7057-0005
	customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO ₂ , O ₂ – or pressure, depending on unit	Hard copy inside folder	7007-0005
Qualification documents	IQ/OQ documents – supporting documents for validation performed by customers,	Digital in PDF format	7057-0001
	consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO ₂ , O ₂ , pressure, depending on unit	Hard copy inside folder	7007-0001
Shelf	Incl. shelf supports	stainless steel	8012-1902
WLAN kit	The kit contains a client bridge that establishes a wireless connection between BINDER units and APT-COM4, LIMS, or customer-specific software via an Ethernet interface. It provides an alternative solution in situations where units need to be positioned in locations without an on-site Ethernet connection. In secured networks, installation and configuration must be performed by the customer's IT service.		8012-2262
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	_	8012-2250



SERVICES

Designation	Description Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts		Article Number DL20-0710	
BRONZE 3-year maintenance contract				
C 1:1 1: 1:C 1	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate		0101	
Calibration certificate, temperature	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate		0102	
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder	DL41-0	0210	
Execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder	DL44-0	0500	
GOLD 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, replacement of wear parts, calibration of one test temperature specified by the user in the center of the usable space, including certificate	DL20-0	0910	
Maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate	DL20-0	0200	
SILVER 3-year maintenance contract	of control response 20% discount on spare parts, festing of all key functions, calibration of one fest		0810	
Temperature measurement, 20 measuring points	Temperature measurement with 20 measuring points with a set value specified by the user, including certificate	DL30-0	0620	
Unit commissioning	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work)	DL10-0	2300	
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming	05 DL10-0	2700	



NOTES

Of Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

BINDER GmbH

Tuttlingen, Germany
TEL +49 7462 2005 0
FAX +49 7462 2005 100
info@binder-world.com
www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd.

Kowloon, Hong Kong, P.R. China TEL +852 39070500 FAX +852 39070507 asia@binder-world.com www.binder-world.com

BINDER Inc.

Bohemia, NY, USA TEL +1 631 224 4340 FAX +1 631 224 4354 usa@binder-world.com www.binder-world.us

BINDER Environmental Testing Equipment (Shanghai) Co., Ltd.

Shanghai, P.R. China TEL +86 21 685 808 25 FAX +86 21 685 808 29 china@binder-world.com www.binder-world.com



Helago-CZ s.r.o. Kladská 1082 500 03 Hradec Králové Tel.: 495 220 229, 495 220 394

Fax: 495 220 154 http://www.helago-cz.cz E-mail: info@helago-cz.cz